
NOOĀNCE

P A R I S

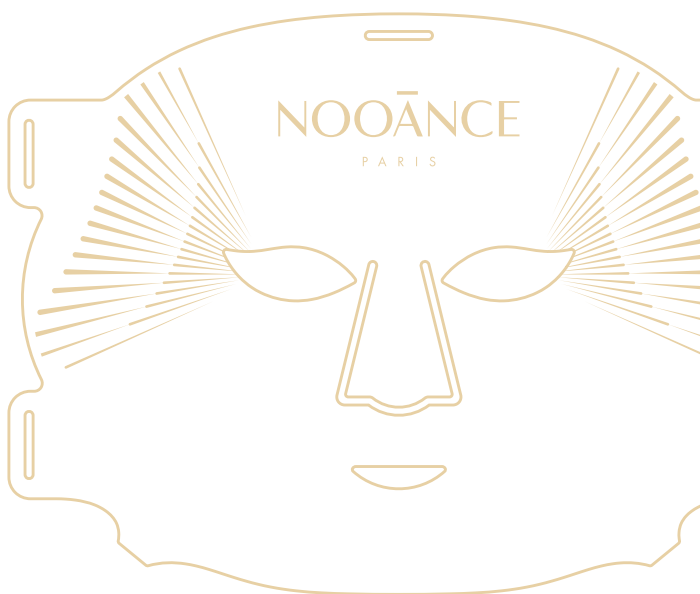MASQUE LED VISAGE
LED FACIAL MASK

PRO

Français ————— 03-19

English ————— 20-37

Deutsch ————— 38-55

Español ————— 56-73

Italiano ————— 74-91

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Sommaire

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Table of Contents

Recommendation of Use	21
Safety instructions	
Cleaning	
Contraindications	
Warnings	
Precautions	
Potential side effects	
Disclaimer	26
Product Representation	26
Customer Service	26
The Science of Light Therapy	27
Product Information	28
Device Specifications	
Contents of the Box	
Product Overview	30
Usage	31
Before First Use	
First Use and Preparation	
User Guide	31
Battery Operation	
Use of Mains Power	
Frequency of use	
Storing	
Warranty	
Air Travel	
Disposal	
Troubleshooting	35
Symbols	36

Recommendation Of Use

Safety instructions

- Read all instructions carefully before using this device.
- All instructions and this manual are available on the website www.nooance-paris.com or upon request to customer service at bonjour@nooance.com
- DO NOT charge or use the device near water or if it is at risk of coming into contact with water. Do not immerse the device.
- DO NOT use the device if it is not working properly, has been dropped, or damaged. For warranty and repairs, contact customer service.
- DO NOT attempt to repair the device. Repairing or opening the device will void your warranty. The device contains no user-serviceable parts. This device contains non-replaceable batteries.
- DO NOT perform maintenance or disassemble the product. In case of malfunction or damage, contact customer service at bonjour@nooance.com
- DO NOT use the device if it is damaged or the cable is damaged. The device must only be powered by a low voltage of 5V DC. If any part of the cable is damaged, stop using it. Contact customer service at bonjour@nooance.com
- DO NOT use or charge the device with a damaged or defective power adapter.
- DO NOT excessively bend or mistreat the USB charging cable. You could damage the device.
- The device must only be powered by a low voltage of 5V DC. If any part of the main cable is damaged, STOP using it. Contact customer service to replace the cable.
- Use only the provided USB cable. Using an unsuitable charging cable may damage the device.

WARNING: To charge the battery 3.7V, 3600mAh, 13.32Wh, use only the removable power unit supplied with your NOOANCE device.

Cleaning

- Clean the device with a soft, dry cloth.
- DO NOT use solvents or abrasive cleaners on the device.
- The device should not be cleaned or charged by children.

The mask is reusable, and it needs to be cleaned after each use. Follow the instructions below to clean your device.

- 1) Before cleaning, please turn off the device and unplug the power.
- 2) Use a clean cotton cloth moistened with water to clean the casing and ensure that no liquid enters the interior of the device. Stubborn stains can be removed by wiping gently with 70% alcohol.
- 3) Wipe the goggles with 70% alcohol wipes after each use.
- 4) Visually check the cleanliness of the device and goggles, make sure there are no dust, impurities, or any other pollutant, otherwise, please repeat the above step again.
- 5) Drying: please dry the device in the air naturally.

Note: Never immerse the device in water or other cleaning liquids to avoid damage to the device.

Contraindications

- Consult your doctor before using this device if you are taking medication (over-the-counter or prescription) that may cause light sensitivity, or if you have suspicious or precancerous lesions.
- DO NOT use the device if you have active cancer or a cancerous lesion.
- DO NOT use the device on areas with phlebitis, thrombophlebitis, or varicose veins.
- DO NOT use the device if you have a known allergy to medical-grade silicone. If in doubt about allergies, contact your doctor before using the device.
- DO NOT use the device if you have a medical history of seizures or epilepsy triggered by light.

Perform a light sensitivity test before using the device on your face:

- Turn on your device.
- Place the device on your arm and use it for 3 minutes.
- After the test, visually inspect the treated area. If the area turns red and the redness persists for more than 2 hours, your skin is sensitive to light. Stop using the device and contact customer service at bonjour@nooance.com

Warnings

- DO NOT allow children under 12 years old to use or play with this device.
- This device can be used by individuals with reduced physical, sensory, or mental abilities, or those lacking experience or knowledge, provided they are properly supervised or have been instructed on the safe use of the device and understand the risks involved.
- If you are pregnant or think you may be, DO NOT use the device without consulting a doctor.
- DO NOT use the device if it is wet. If the device has come into contact with water, do not touch it without unplugging it first.
- STOP using the device and consult your doctor if you experience any adverse reactions.
- STOP using the device if you experience symptoms of low blood pressure (dizziness or lightheadedness) during use.
- STOP using the device if it overheats. Turn off the device immediately.
- STOP treatment if the device feels hot or if your skin emits an uncomfortable warmth.
- DO NOT let children play with the device. Cleaning and charging the device should not be done by children.
- DO NOT use the device if you experience headaches caused by light.
- DO NOT use the device to target rosacea, moles, warts, open wounds, or any other skin condition without medical advice.
- DO NOT use the device on depigmented spots or areas of skin affected by vitiligo.

- DO NOT point the LED light directly into your eyes or another person's eyes. You may experience temporary visual impairment.
- DO NOT use the device if you have open wounds.
- DO NOT use the device with unauthorized accessories.
- DO NOT use solvents or abrasive cleaners on the device.
- DO NOT use the device for purposes other than those specified in this manual or with accessories not authorized by the manufacturer.
- DO NOT use the device if your skin is sunburned. Wait until your skin heals before using it.
- DO NOT use the device if you are receiving medical treatment for conditions involving photosensitivity.
- DO NOT use the device if you are on medication that causes photosensitivity.
- DO NOT use the device if you have had recent surgery or other surgical procedures.
- DO NOT use the device on bacterial infections as this could worsen the symptoms.
- DO NOT use the device on painful areas. Consult your doctor if you experience pain during use.
- DO NOT use the device on individuals under 18 years of age without medical advice.
- No modifications to this equipment are allowed.
- Contact your doctor if you have any medical questions regarding the use of this device.

Precautions

- DO NOT scratch or damage the light-emitting surface.
- Your warranty will be void if you do not use the device in accordance with the instructions in the manual.
- DO NOT use the device outdoors to avoid damage.
- DO NOT place the device in direct sunlight.

- DO NOT sit, lie down, or step on the device to avoid breaking it.
- DO NOT use or store the device near hot surfaces.
- Store the device according to the protocol described in this manual.
- The device is designed for cosmetic use only by healthy individuals. Please contact your doctor if you are unsure of your health condition.
- The device's battery is not removable by the user. If the device needs to be disposed of, do not burn it. Dispose of the device in accordance with local regulations.
- Inspect the device and do not use it if you see signs of damage, such as battery leaks. Contact our customer service.
- Consult your doctor before using the device if you have a fever, acute inflammation, severe circulatory problems, light sensitivity, anesthesia, or desensitized areas.
- We strongly recommend exercising caution when taking medications. If in doubt, consult your doctor.

Note: For hygiene reasons, we advise not sharing your device with others without disinfecting the part in contact with the skin using 70% alcohol.

Potential side effects

- Immediately after your session, the skin may appear red. The skin should return to normal within 1 to 2 hours.
- Temporary headaches may occur after light therapy. Consult your doctor if headaches persist for an extended period.

Disclaimer

- This product is not intended to diagnose, treat, or prevent any disease. NOOANCE makes no claims or warranties regarding this product's ability to treat any mental, physical, or dermatological condition. Consult a qualified healthcare professional for any condition requiring medical attention.
- The NOOANCE Face Mask is not intended for use in the diagnosis of diseases and other conditions, nor for the cure, mitigation, treatment, or prevention of diseases, and is not intended to alter the structure or function of the body.

Product Representation

NOOANCE continually improves its products using the latest technologies. As a result, the actual product may differ slightly from the product shown in the images on the packaging and in this user guide.

Customer Service

For more information on using or maintaining the device, or to obtain additional information, please contact Customer Service via email at bonjour@nooance.com.

The Science Of Light Therapy

The skin has the remarkable ability to absorb light energy. It uses this energy to boost the cells to produce more energy (ATP), which increases their performance and allows them to carry out their natural functions, but at a faster rate.

Cells supercharged by light therapy use this energy to boost the skin's fibroblasts to produce more collagen and elastin proteins and to stimulate the natural processes of cell renew.

The benefits of light therapy :

- boost cell regeneration
- stimulate skin metabolism
- repair and soothe irritations

The NOOANCE Professional LED Mask features 198 LEDs with a precision of ± 2 nanometers :

- 99 red LEDs with a wavelength of 633 nanometers
- 99 near-infrared LEDs with a wavelength of 830 nanometers.

Red and near-infrared LEDs help to :

- reduce wrinkles
- improve firmness
- give an immediate radiance boost


A powerful synergy of 633 nm and 830 nm wavelengths delivering visible results to revitalize cells and improve your skin's radiance, firmness and health.

Product Information



Device Specifications

Product name	LED FACIAL MASK PRO
Model	MK-99A
LED wavelengths	633 nm ± 2 nm (Red) 830 nm ± 2 nm (Infrared)
Max. output optical power density	50mW/cm²
Energy delivered	30J/cm²
Lamp Classification	Exempt group (Risk group according IEC 62471: 2006)
Operating temperature	0°C to 40°C 32°F to 104°F
Storing and transport temperature	-5°C to 55°C 23°F to 131°F
Operating humidity	30% to 95% without condensation
Storage and Transport humidity	20% to 95% without condensation
Charging supply	5V=2A
Charging port	USB-C
Charging cable	USB-C to USB-A
Charging cable length	2 meters / 6.562 feet
Charging adapter power output	2A
Battery	3.7V, 3600mAh
Battery type	Lithium-ion
Charging time	240 minutes / 4 hours
Recommended use	30 to 60 minutes per week
Timer	10 minutes
Dimensions of the device	290*216*9 mm
Weight of the device	0.512kg / 18.06oz

Contents Of The Box

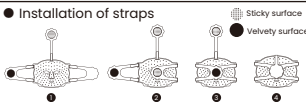Name	Package Contents	Additional information
LED FACIAL MASK PRO		● Clean the device only with a clean, damp cloth or water wipe. You can also use an antibacterial cleaner or wipe. ● DO NOT clean the appliance when it is switched on.
Adjustable Straps		● Installation of straps 
Power Supply		● Only use the provided equipment and avoid fast chargers as they may damage the battery.
Adapter		● Adapter EU / US / UK ● Please use the adapter to connect the controller for charging.
USB Cable		● Please note that the USB cable is only for connecting the LED mask, not for charging the controller.
Removable eye protection		● When using an LED mask, red light and infrared light pose no danger to the eyes. ● The red light of the LEDs may appear orange after a few seconds of use as the eye adapts to this brightness. This is completely normal. ● It is important not to use silicone mask protectors and silicone eye protectors at the same time. Using them together can cause uncomfortable pressure on the eyes.
Silicone Eye Protectors		● The use of goggles is recommended in the event of discomfort or hypersensitivity to light. ● If the silicone eye protectors are adequate, they are sufficient; if the light remains too strong, it is preferable to opt for the dark glasses alone under the mask.
Controller		● Fully charge the controller for 4 hours before using for the first time. To maximise the life of the lithium battery, we recommend that you fully charge and discharge it during normal use. ● The indicator light flashes during charging ● Unplug the controller from the mains before using it. The device does not work when it is charging. ● Before charging the controller, make sure that the plug and socket are completely dry. Failure to do so may result in an electric shock, short-circuit or fire.
Storage Bag		● The device should be stored at room temperature in a dry place and in its original packaging to prevent it from deforming, stretching or being compressed. Avoid dusty, sunny or damp locations.
Manual		● Read this manual carefully before using your mask.
Waterproof Plug		● For safety reasons, please use the waterproof Plug to block the charging port when not charging.

Product Overview

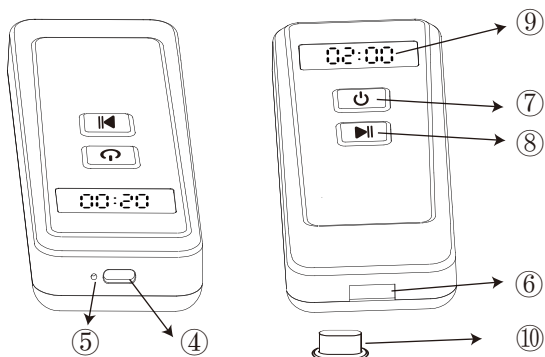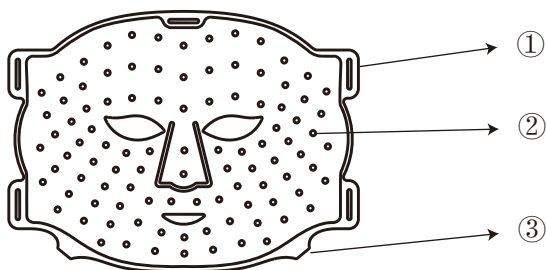- ①.Adjustable Straps Connection Port
- ②.Led -Emitting Surface
- ③.USB Cable Connection Port
- ④.Charging Port
- ⑤.Charging & Operating Indicator Light
- ⑥.USB Cable Type-C Connection Port
- ⑦.Power ON/OFF Button
- ⑧.Pause/Start Button
- ⑨.Operating Timer Digital Display
- ⑩.Waterproof Plug

Usage

Before First Use :

We strongly recommend that users carry out a light sensitivity test before first use.

1. Place the device on the inside of the forearm.
2. Switch on the device to start your session.
3. Turn the mask off after 5 minutes (the timer is set for 10 minutes by default).

If your skin becomes excessively hot or itchy during the session, stop immediately. You may be sensitive to light.

Check your skin regularly for the next 24 hours. If your skin is red and/or itchy, the appliance is not suitable for you. Please contact customer service.

First Use and Preparation

1)Unpack the package and check off all the parts described in section "Contents of the box".

2)Before use it is important to fully charge the controller.

-The USB-C port of the USB cable is plugged into the adapter, and the other port(type-c) is connected to the Type-C Connection Port of the controller. Plug the adapter into the power outlet (100-240Va.c., 50Hz/60Hz) to power supply the device. When charging, the battery indicator will flash white. When fully charged, the indicator lights will stop flashing and change to a constant white light


3)When the battery level is less than 10%, the indicator lights flash white. To recharge the device, follow the instruction above.

Note:


- Once fully charged, detach the controller from the charging.
- Use the adapter (input: 100-240Va.c., 50Hz/60Hz, output: DC 5V, 2A) supplied with the unit. The AC adapter used should be certified according to IEC/EN 61558 or IEC/EN 62368-1.

User Guide

Preparation: Before starting the treatment, please thoroughly cleanse your skin.

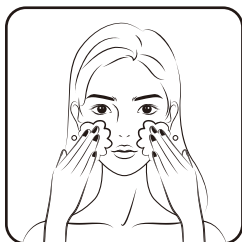
Insert the Removable eye protection: Align the narrow edge of the eye protection with the eye socket of the mask's inside, and push it into the outside of the mask. The broader edge of the eye protection should be facing the user's skin.

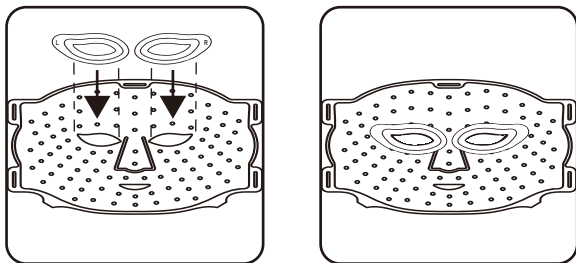Note: If the user wants better eye protection, they can choose to use Silicone Eye Protectors.


Connect The Adjustable Straps

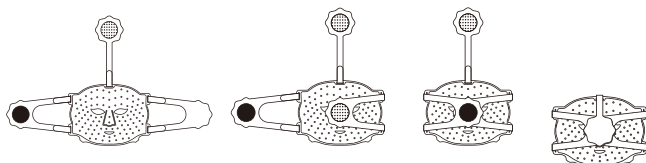Place the mask over your face: Apply the mask to your face, and adjust the straps to ensure that the mask sits well on your face / desired position.


Tip: To optimise the results of your Nooance mask, we recommend that you use one of our aqueous serums recommended on our website www.nooance-paris.com. At the end of the session, you can also apply a moisturising cream or serum.

Connect The Mask To Controller

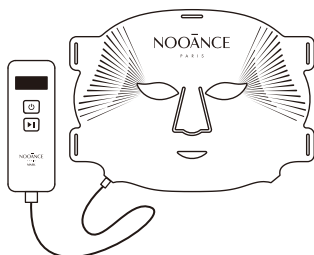
Connect the mask and the controller with the USB-cable provided.

Battery Operation

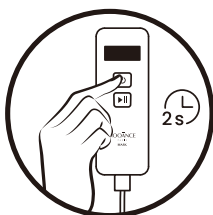Press and hold the ON/OFF button for 2 seconds to turn on the mask on the controller.

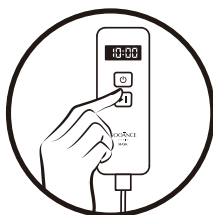Press the "Start/Pause" button, the indicator light will turn on and start the 10-minutes session. Once the session is complete, the mask will turn off automatically.

Use Of Mains Power

Connect the mask directly to the mains using the 2m cable supplied and the mains adapter (you will need to time your session manually).

Do not connect the controller/battery in this case.

Frequency Of Use

We recommend using the mask for 30 minutes once a week (or 3 x 10 minutes spread over the week) for results on fine lines and wrinkles from as early as 16 weeks. In the first month, you can start by doubling the time per week if you want to speed up the cell renewal process and boost skin firmness.

Storing

- 1) Avoid direct sunlight on the LED FACIAL MASK PRO while storing.
- 2) It is recommended to store the LED FACIAL MASK PRO in its original packaging.
- 3) Store in an environment with a temperature of $-5\sim 55^{\circ}\text{C}$, a humidity of 20-95%RH and an atmospheric pressure of 700hPa-1060hPa.

Warranty

- The LED mask comes with a 24-month guarantee for the original purchaser, subject to proof of purchase from an authorised dealer. This guarantee covers normal use with approved accessories. In the event of a fault, contact our customer service department.
- Please note that the guarantee does not cover damage caused by accidents, incorrect use, use of non-recommended accessories, or poor electrical connection.

Air Travel

- The device is allowed in the cabin. If you wish to check it into the hold, take the controller/Battery with you in the cabin. The main airlines refuse to accept lithium batteries in the hold. If in doubt, consult your airline.

Disposal

DO NOT dispose of the device with household waste. The device should be recycled. The device contains a Lithium-ion battery. For protection of the environment, please safely dispose of the device following your local guidelines.

Troubleshooting

If the actions suggested fail to fix the problem, or the problem reoccurs, and your device is within the warranty period please contact our customer service by writing to bonjour@nooance-paris.com.

⚠ DO NOT attempt to investigate any further.

⚠ DO NOT attempt to modify the device in any way. The device has NO parts you can repair yourself.

Problem	Cause	Solution
The device doesn't switch on.	The device is not fully charged, because the charging cable was not properly inserted into the power outlet.	Insert the charging cable properly into the power outlet.
	The device is not fully charged, because the charging cable was not properly inserted into the device.	Insert the charging cable properly into the device and charge the device.
	The power outlet didn't work during charging.	Try a different 5V adapter or power outlet.
	The on/off button has not been pushed 2 seconds full second.	Push the on/off button for 2 seconds straight.
The device switches off unexpectedly.	The battery level is too low.	Charge the Controller.

Symbols

Symbol	Explanation
	Indicates the device is CE Marked and authorized to sell in Europe. (Lower-Voltage-Directive)
	The Waste of Electrical and Electronic Equipment (WEEE), Indicates the device may not be discarded as household waste (2012/19/EU).
	Indicates that the electromagnetic radiation from the device is below the limits specified by the Federal Communications Commission and the manufacturer has followed the requirements of the Supplier's Declaration of Conformity authorization procedures.
	This warning symbol on the device label is a warning to read the enclosed documents before use and take all safety precautions required.
	Consult instructions for use.
	Date of Manufacture.
	This symbol means keep the device dry.
	This symbol indicates the serial number of the device.
	This symbol refers to the batch code of the device.
	This symbol indicates the USB connector.
	This symbol means that the device uses a Lithium-ion Polymer battery.
	This symbol means keep the device away from direct sunlight.
	This symbol means the upper and lower operating temperatures that the device can be safely used at.
	This symbol indicates the upper and lower operating humidity that the device can be safely used at.

	This symbol means that the cardboard packaging that the device came in can be recycled.
	IP22
	Unique Device Identifier.
	This symbol implies conformity to UK Directives on health, safety and environmental protection standards. UKCA stands for UK Conformity Assessed.
	Indicates the device is CA65 Marked and authorized to sell in the United States.
	Indicates the device is ROHS Marked and authorized to sell in Europe.
	Class III equipment
	To indicate that a specific separate supply unit is required for connecting electrical equipment to the supply mains.
	Manufacturer
	Country of manufacture
	Importer
	EU Representative

Shenzhen Oren Medical Technology Co., Ltd.
Address: 101, Building 3, No. 40, Fuxin Street, Huaide Community,
Fuyong Street, Baoan District, Shenzhen, China

Country of manufacture
Made in China

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