
NOOANCE

PARIS

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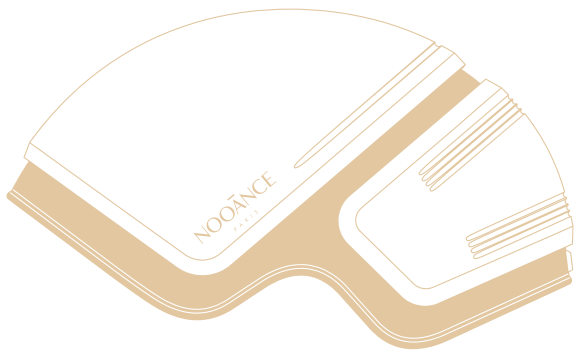HAIR HELMET

DOUBLE TECHNOLOGY LASER & LED

M-282 PRO

Français ————— **03-26**

English ————— **28-51**

Sommaire

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NOOÅNCE LASER AND LED HELMET
 Model: M-282 PRO
 Input: 5V = 2A
 Applied Part: Class 3R Laser Product
 Use only the specified adapters and USB cables.

NOOÅNCE INNOVATIONS
 22 Rue Beaujon, 75008
 Paris, France

221115

LOT 4722OR

SN M-28200001

(01)6970006970918
 (11)221115
 (10)4722OR
 (21)M-28200001

MADE IN CHINA

LASER RADIATION
 AVOID DIRECT EYE EXPOSURE
 CLASS 3R LASER PRODUCT
 Wavelength: 650±5nm,
 Output: $+5mW, CW$

IEC 60825-1:2014 SAFETY OF LASER PRODUCT

NOOÅNCE LASER AND LED HELMET
 Model: M-282 PRO
 Input: 5V = 2A
 Applied Part: Class 3R Laser Product
 Use only the specified adapters and USB cables.

NOOÅNCE INNOVATIONS
 22 Rue Beaujon, 75008
 Paris, France

221115

LOT 4722OR

SN M-28200001




(01)6970006970918
 (11)221115
 (10)4722OR
 (21)M-28200001

MADE IN CHINA

LASER RADIATION
 AVOID DIRECT EYE EXPOSURE
 CLASS 3R LASER PRODUCT
 Wavelength: 650±5nm,
 Output: $+5mW, CW$

IEC 60825-1:2014 SAFETY OF LASER PRODUCT

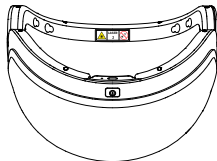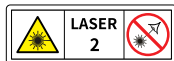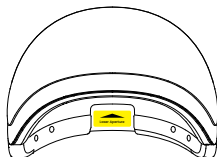NOOÅNCE INNOVATIONS
 22 Rue Beaujon, 75008 Paris, France

Conçu en France
 Fabriqué en Chine

Table of Contents

Contraindications	29	Charging	42
Precaution	29	Care And Storage	43
Warnings	30	Cleaning	43
Device Specifications	33	Troubleshooting	44
Device Accessories		Electrical Safety	46
Light Sensitivity Test	35	Travel and International Use	47
How To Use The Device (Two ways to use)	36	Disposal	47
General Use		Warranty	48
Visibility of continuous mode and pulse mode		Symbols	49
Installation of Sponge Pads			
Session Protocol			

Contraindications

- Do not use the NOOANCE helmet if you have a medical condition that makes your skin sensitive to sun or light or if you are taking photosensitizing medications.
- Do not use the NOOANCE helmet if you have cancer, skin diseases, or other conditions on your scalp.
- The NOOANCE helmet is only for people over the age of 18.
- Please consult your doctor prior to use if you are, or suspect you are, pregnant or breastfeeding.

Precaution

- If you are taking any medications or products that make you sensitive to light, you should test your NOOANCE helmet to make sure you do not have a reaction. Some of the things that can cause increased sensitivity to your NOOANCE helmet could be medicines for colds, allergies, pain medicine and drugs to treat infections.
- ⚠ Only use your NOOANCE helmet as directed. Do not try to change the way your NOOANCE helmet operates. Altering controls, adjustments, performance, or procedures other than those specified herein may result in exposure to hazardous radiation (laser). Never use acetone or other solvents on any part of your NOOANCE helmet. The use of solvents on the window can block the laser light from coming through.

Warnings

- Before using this NOOĀNCE helmet, you must take a skin sensitivity test.
- Use the NOOĀNCE helmet with caution, as handling in a manner other than that specified in the NOOĀNCE manual may result in dangerous radiation exposure.
- Do not stare directly at the red laser light or its reflection in a mirror because it can temporarily irritate your eyes. If you stare at the light for too long, it could harm your eyes. Also, never use a magnifying glass when using the NOOĀNCE helmet because it can cause temporary irritation to your eyes.
- The NOOĀNCE helmet contains Class 3R lasers that may cause harm to the retina and eye irritation if the lasers are pointed directly at the eyes. 
- Stop using NOOĀNCE helmet if you experience itching or tingling of the scalp. If the problem persists one hour after the end of the session, stop use and consult your doctor.
- Keep the NOOĀNCE helmet out of reach of children and pets.
- Do not place the NOOĀNCE helmet in a wet location.
- Do not wet the helmet, controller (battery), plug or adapter. This is an electrical hazard that can result in personal injury. 
- Please do not use the NOOĀNCE helmet while charging the controller (battery). This can lead to faults or electrical hazards. 


- Do not pull or forcibly detach the power cord when using the NOOANCE helmet. This could cause a fault or an electrical hazard.
- Please do not use it in the bathroom or the shower. Please do not use a NOOANCE helmet that has been wet with water. 
- Keep the NOOANCE helmet away from hot objects like ovens and indoor heaters.
- Make sure your hair is dry when using the device. Do not use the NOOANCE helmet when your hair is wet.
- Use only the controller (battery) and plugs/adaptors provided with the NOOANCE helmet.
- Use the NOOANCE helmet according to the instruction manual provided. Please do not use it for any other purpose.
- If you drop the NOOANCE helmet or suspect that the NOOANCE helmet is faulty (e.g., the LED or laser does not work, the NOOANCE helmet is cracked). Please get in touch with customer service for repair or replacement.
- Do not use the NOOANCE helmet when room temperature exceeds 33.3°C(92°F). This could cause dryness, redness or irritation of your scalp.
- If the cord is damaged on your NOOANCE helmet, do not use the NOOANCE helmet because you could get an electrical shock. Please get in touch with customer service for repair or replacement if the cord gets damaged.

- Do not attempt to modify or repair the NOOANCE helmet.
- Always unplug the NOOANCE helmet from the power outlet prior to cleaning the NOOANCE helmet. Failure to do so may result in an electrical hazard that may cause harm or injury.


Device Specifications

Product Name	LASER AND LED HELMET
Model No.	M-282 PRO
Output optical power density	<5mW
Wavelength	650±5nm
Number of Light Emitting Diodes (LED)	200
Number of Lasers (Class 3R)	82
Mode	Mode 1 is continuous mode; Mode 2 is pulse mode.
Operating temperature/humidity/pressure	5°C-35°C, 15%-90% RH, 700 hpa to 1060 hpa
Storage & transportation temperature/humidity	-25°C-70°C, 15%-90% RH
Ingress of water or particulate matter	IP22 (Enclose equipment without protection against ingress of water)
Charging supply	5V = 2A
Battery	3.7V, 3600mAh (Lithium-ion)
Session time	25 minutes
Charging time	4 hours
Recommended use	Three times a week on alternate days
Dimensions of the device	259.94*218.12*155.86 mm
Weight of the device	0.585kg / 20.63oz


Device Accessories

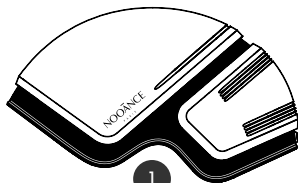1

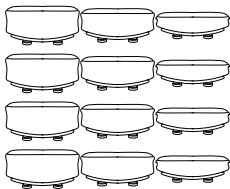2

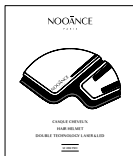3

4

5

6

7

8

9

10


No.	Accessories	Quantity
1	LASER AND LED HELMET (M-282 PRO)	1
2	Sponge Pad (Size: 1.5cm, 2.5cm, 3cm)	12
3	User manual	1
4	Controller (battery)	1
5	Power Supply	1
6	EU Adapter	1
7	US Adapter	1
8	UK Adapter	1
9	USB Cable	1
10	Storage Bag	1

Light Sensitivity Test

We highly recommend users perform a light sensitivity test before their first use of the NOOANCE helmet.

1. Place the NOOANCE helmet on the inner forearm.
2. Turn on the NOOANCE helmet to start a session.
3. Wait until the session is done. The session will take 5 minutes to complete.
4. If your skin feels hot or uncomfortable during the session stop immediately! You might be sensitive to light.
5. The NOOANCE helmet will automatically turn off after the session is done. Remove the NOOANCE helmet and wait for at least 24 hours.

Check your skin regularly over those 24 hours. If your skin appears red or itchy, the NOOANCE helmet is unsuitable for you. In the case of light sensitivity, please get in touch with customer support. If you do not see a reaction, proceed with the session instructions.

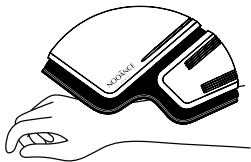How To Use The Device


(Two ways to use)

General use

Preparation: Make sure your hair is dry before using the NOOANCE helmet.

① Connecting to the controller (battery)

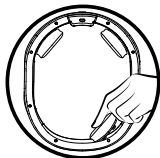Step 1

Please select the appropriate sponge pad size to fit inside the NOOANCE helmet.


Step 2

Connect the NOOANCE helmet to the controller (battery) using a USB cable.

Step 3


Position the NOOANCE helmet in its upright position on the top of your head.

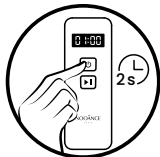Step 4

Press and hold the ON/OFF  button for 2 seconds to switch the device on.

(**Note:** If the start/pause button is not pressed within five minutes, the NOOANCE helmet will turn off automatically.)

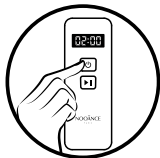Step 5

Gently press the ON/OFF  button to select your preferred mode. (Mode 1 is continuous; Mode 2 is pulsed.)

Step 6


Press the start/pause button to begin the 25-minute session. A "beep" indicates the session has started. Once the session is complete, the NOOANCE helmet will turn off automatically with two "beep" sounds.

IMPORTANT: Unless the NOOANCE helmet is positioned on the user's head, it will not begin a session. If the NOOANCE helmet is removed from the head during a session, it will pause automatically. To resume the session, position the NOOANCE helmet on the user's head within 5 minutes and press the start/pause button.

Step 7

If you need to pause the session at any time, press the start/pause button, and the timer will stop. Press the start/pause button to resume the session. A start/pause audible "beep" sound will indicate the session has resumed. If the start/pause button is not pressed within 5 minutes during a paused session, the NOOANCE helmet will turn off automatically.

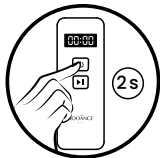Step 8

To turn off the NOOANCE helmet, hold the button for 2 seconds.

IMPORTANT: The NOOANCE controller (battery) has a memory in case of a power failure. In the event of a power failure or accidental disconnection of the USB cable during use, simply press the start/pause button again. The helmet will resume the timer from where it left off. To switch off the device, press the ON/OFF  button for one second. Then follow steps 2 to 8.

② Connect to the mains plug with adapter

Step 1

Connect the helmet to the mains using the correct adapter.


Step 2

Press and hold the On/Off  button to turn on the NOOANCE helmet default start mode 1, and the indicator light on the NOOANCE helmet will always be white.

Step 3

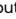Press the On/Off  button again to switch to mode 2, and the indicator light on the NOOANCE helmet will flash white.

Step 4

Place the NOOANCE helmet vertically on your head after selecting your preferred mode. Once the session is complete, the NOOANCE helmet will turn off automatically with two "beep" sounds.

Step 5

To turn off the NOOANCE helmet, hold the button for 2 seconds.

Visibility of continuous mode and pulse mode

- Figure 1: Continuous mode. The LED and laser light up continuously, and the oscilloscope will display a fixed frequency.
- Figure 2: The pulse mode frequency is set to 1.97kHz, so the rapid flicker of the laser is invisible to the naked eye. But it shows up clearly on the oscilloscope.

Figure 1: Mode 1 is continuous mode

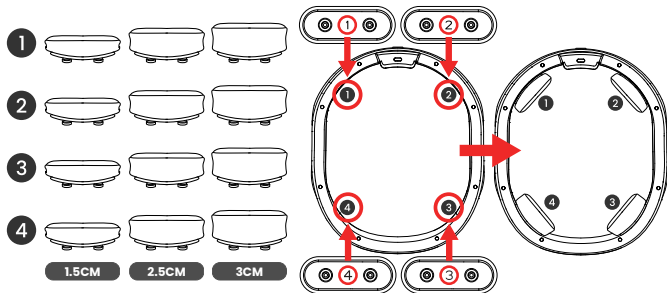
Figure 2: Mode 2 is pulse mode

Invisible pulse mode (2) is used in most laser machines or LEDs dedicated to professionals. It has the disadvantage of producing background noise associated with the pulse, which depends on each person's auditory sensitivity, which is why we have also created a continuous mode so that everyone can use our helmet.

Installation of Sponge Pads

There are 3 sizes to choose from. The corresponding number is located below each pad. Please put the correct number in the specified position.

Session Protocol

You should use the NOOANCE helmet for 25 minutes, every other day, on non-consecutive days.

Use the NOOANCE helmet for 16 weeks. Once the 16-week period is completed, use as frequently as needed to maintain the desired results without exceeding every other day for 25 minutes. Discontinuing the use of the NOOANCE helmet may result in the resumption of hair loss.

Charging

Fully charge your NOOANCE helmet for 4 hours before your first use. Only use a 5V adapter or computer USB port (avoid fast chargers as they can damage the battery). To charge the NOOANCE helmet insert the USB-C cable into the controller (battery). Insert the USB-A side of the cable into a USB power supply. When the controller (battery) is charged, the LED indicator on the device flashes white, and when the controller is fully charged, the indicator will be steady white. Before charging, make sure that the plug and socket are completely dry. Failure to do so may result in electric shock, short circuit, or fire.

Note: Always unplug the NOOANCE helmet before use, the NOOANCE helmet will not work while being charged.

Care And Storage

- Make sure that the NOOANCE helmet is not bent, stretched or crushed during storage or transport.
- Keep the NOOANCE helmet away from sunlight and store the NOOANCE helmet in a dry location (e.g. away from windows, sinks, bathtubs, kitchens, and bathrooms) at room temperature after each use.

Cleaning

- The NOOANCE helmet should be unplugged prior to any cleaning.
- Use a dry microfiber cleaning wipe or cloth to clean the device.
- For stains that cannot be removed with a dry cloth, a lightly water-moistened cloth may be used to clean the NOOANCE helmet.
- Do not use solvents to clean the NOOANCE helmet.

Troubleshooting

If the actions suggested fail to fix the problem, or the problem reoccurs, and your device is within the warranty period please return to NOOANCE.

⚠ DO NOT attempt to investigate any further.

⚠ DO NOT attempt to modify the NOOANCE helmet in any way. The NOOANCE helmet has NO parts you can repair yourself.

Problem	Cause	Solution
The NOOANCE helmet doesn't switch on.	The controller (battery) is not fully charged because the charging cable was not properly inserted into the power outlet.	Insert the charging cable properly into the power outlet.
	The remote control is not fully charged because the charging cable was not properly inserted into the power outlet.	Insert the charging cable properly into the NOOANCE helmet and charge the NOOANCE helmet.
	The power outlet didn't work during charging.	Try a different 5V adapter or power outlet.
	The user does not wear the NOOANCE helmet upright on his/her head.	Unless the NOOANCE helmet is positioned on the user's head, it will not begin a session.
	The NOOANCE helmet connects not only to the controller but also to a power outlet.	After the controller (battery) is powered on, it is in charge mode by default. When the NOOANCE helmet cannot be activated in charging mode, disconnect the power, and the controller (battery) will resume regular use.

	Allow NOOANCE helmet internal software to reset itself.	Hold for 3-5 seconds and allow the NOOANCE helmet's internal software to reset itself.
The NOOANCE helmet is powered on but pressing does not activate the NOOANCE helmet.	The NOOANCE helmet will not power on unless the electronic sensor detects that the NOOANCE helmet is positioned on the users head. This is a built-in safety feature to protect the eyes from direct laser light exposure.	Activate the NOOANCE helmet by placing it vertically over your head.
The NOOANCE helmet is not lighting up, but the controller is on.	The power of the controller is low.	Please fully charge the controller before use.
	If the device still fails to light up after following the instructions and troubleshooting.	Please get in touch with customer service.
The NOOANCE helmet controller doesn't show the numbers clearly on the timer.	Try shutting down the device, restarting it, and allowing the NOOANCE helmet's internal software to reset itself.	If the timer cannot be displayed clearly after the NOOANCE helmet is restarted, please get in touch with customer service.

Electrical Safety

To reduce the risk of hazard by electric shock:

- Keep away from water. Contact with water may result in an electrical shock.
 - Do not use while bathing.
 - Do not place or store your NOOANCE helmet where it can fall or be pulled into a tub or sink.
 - Do not expose your NOOANCE helmet to running water, or drop into water or any other liquid. If your NOOANCE helmet falls into water, unplug it immediately from the power outlet.
 - If using your NOOANCE helmet near water, the use of an electrical outlet protected by a Ground Fault Circuit Interrupter is recommended.
 - Do not dispose of your NOOANCE helmet in a fire or store in high temperatures greater than 140°F (60°C).
 - Do not turn on or use your NOOANCE helmet while it is charging.
 - The power supply cannot be separated during use. Do not position the equipment in a manner that will make it difficult to plug in or remove the power cord.
 - The power cord is very long. To avoid danger, do not wrap the cord around your neck.
- ⚠ Basic safety precautions should be observed when using this, or any other electrical device, especially when children are present.

Travel and International Use

The NOOĀNCE helmet's lightweight design makes it a perfect travel companion. Here are some travel tips for you:

- If you are travelling to a different region from where you purchased the NOOĀNCE helmet, make sure that the NOOĀNCE helmet is only supplied with a low voltage of 5V DC.
- Confirm with the steward that the NOOĀNCE helmet can be carried and/or used on the airplane prior to flight.

Disposal

DO NOT dispose of the device with household waste. The device should be recycled. The device contains a Lithium-ion battery. For protection of the environment, please safely dispose of the device following your local guidelines.


Warranty


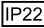
The NOOANCE helmet is covered by a warranty against defects in material and against defects in workmanship and materials for 24 months from the date of purchase under normal conditions of use. This warranty extends only to the original retail purchaser with original proof of purchase and only when purchased from an authorized retailer or reseller. This warranty only extends if the NOOANCE helmet is used in conjunction with authorized Device accessories. If the NOOANCE helmet should become defective within the warranty period, contact the retailer that you purchased the NOOANCE helmet from. This warranty does not cover the NOOANCE helmet damaged by:

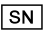- Accident, misuse, abuse, or alteration.
- Use with unauthorized accessories or use other than as instructed.
- Connecting to incorrect current and voltage.

Symbols

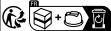Electromagnetic compatibility complies with IEC 60601-1-2 Standard.

symbol	Explanation
	The BF symbol, which indicates this device is in accordance to the degree of protecting against electric shock for type BF equipment.
	Indicates the Device is CE Marked and authorized to sell in Europe. (Lower-Voltage-Directive)
	The Waste of Electrical and Electronic Equipment (WEEE), Indicates the Device may not be discarded as household waste (2012/19/EU).
	Indicates that the electromagnetic radiation from the device is below the limits specified by the Federal Communications Commission and the manufacturer has followed the requirements of the Supplier's Declaration of Conformity authorization procedures.
	This warning symbol on the device label is a warning to read the enclosed documents before use and take all safety precautions required.
	Consult instructions for use.
	IP22
	Unique Device Identifier.
	This symbol means keep the Device dry.

	This symbol indicates the serial number of the device.
	This symbol refers to the batch code of the device.
	This symbol means keep the device away from direct sunlight.
	This symbol indicates the upper and lower operating temperatures that the device can be safely used at.
	This symbol means the upper and lower operating humidity that the device can be safely used at.
	This symbol means that the cardboard packaging that the device came in can be recycled.
	This symbol implies conformity to UK Directives on health, safety and environmental protection standards. UKCA stands for UK Conformity Assessed.
	Indicates the Device is CA65 Marked and authorized to be sold in the United States.
	Indicates the Device is ROHS Marked and authorized to sell in Europe.
	This symbol means that the device uses a Lithium-ion Polymer battery.
	The date of manufacture.

- Labels on package and device are in compliance with IEC 60825-1 Standard:

NOOÃNCE LASER AND LED HELMET
 Model: M-282 PRO
 Input: 5V = 2A
 Applied Part: Class 3R Laser Product
 Use only the specified adapters and USB cables.

NOOÃNCE INNOVATIONS
 22 Rue Beaujon, 75008
 Paris, France


221115

4722OR

M-28200001

(01)6970006970918
 (11)221115
 (10)4722OR
 (21)M-28200001

MADE IN CHINA

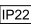LASER RADIATION
 AVOID DIRECT EYE EXPOSURE
 CLASS 3R LASER PRODUCT
 Wavelength: 650±5nm,
 Output: $+5mW, CW$

ICE 60825-1:2014 SAFETY OF LASER PRODUCT

NOOÃNCE LASER AND LED HELMET
 Model: M-282 PRO
 Input: 5V = 2A
 Applied Part: Class 3R Laser Product
 Use only the specified adapters and USB cables.

NOOÃNCE INNOVATIONS
 22 Rue Beaujon, 75008
 Paris, France

221115

4722OR

M-28200001




(01)6970006970918
 (11)221115
 (10)4722OR
 (21)M-28200001

MADE IN CHINA

LASER RADIATION
 AVOID DIRECT EYE EXPOSURE
 CLASS 3R LASER PRODUCT
 Wavelength: 650±5nm,
 Output: $+5mW, CW$

ICE 60825-1:2014 SAFETY OF LASER PRODUCT

NOOÃNCE INNOVATIONS
 22 Rue Beaujon, 75008 Paris, France

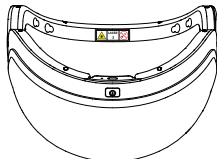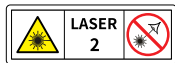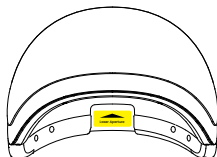Designed in France
 Made in China