
NOOANCE
PARIS

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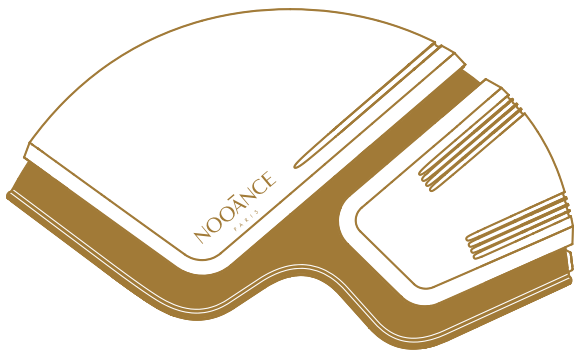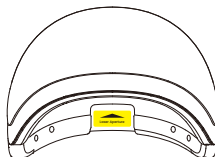
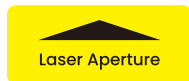
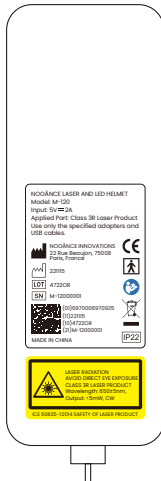
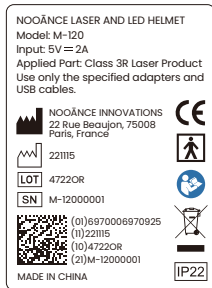
CONÇU EN FRANCE
FABRIQUÉ EN CHINE

Table of Contents

Company Information	36	Troubleshooting	55
Contraindications	37	Electrical Safety	57
Warnings	38	Travel And International Use	58
Precaution	41	Disposal	59
Product Description	42	Warranty	60
Device Specifications	43	Symbols	61
Device Accessories			
Light Sensitivity Test	45		
How To Use The Device	46		
General Use			
Installation of bearings			
Treatment Protocol			
Charging	52		
Care And Storage	53		
Cleaning	54		

Company Information

Thank you for choosing NOOÄNCE LASER AND LED HELMET. It is important that you read and understand the content of this User Manual before using your NOOÄNCE LASER AND LED HELMET.

Company

NOOÄNCE INNOVATIONS

Address

22 Rue Beaujon, 75008 Paris, France

Contraindications

- Do not use the NOOANCE helmet if you have a medical condition that makes your skin sensitive to sun or light or if you are taking photosensitizing medications.
- Do not use the NOOANCE helmet if you have cancer, skin diseases, or other conditions on your scalp.
- The NOOANCE helmet is only for users over 18.
- Please consult your doctor prior to use if you are, or suspect you are, pregnant or breastfeeding.

Warnings

- Before using this NOOĀNCE helmet, you must take a skin sensitivity test.
- Use the NOOĀNCE helmet with caution, as handling in a manner other than that specified in the NOOĀNCE manual may result in dangerous radiation exposure.
- Do not stare directly at the red laser light or its reflection in a mirror because it can temporarily irritate your eyes. If you stare at the light for too long, it could harm your eyes. Also, never use a magnifying glass when using the NOOĀNCE helmet because it can cause temporary irritation to your eyes.
- The NOOĀNCE helmet contains Class 3R lasers that may cause harm to the retina and eye irritation if the lasers are pointed directly at the eyes. 
- Stop using the NOOĀNCE helmet if you experience itchy or tingling scalp. If the problem persists one hour after the end of treatment, discontinue use and consult your doctor.
- Keep the NOOĀNCE helmet out of reach of children and pets.
- Do not place the NOOĀNCE helmet in a wet location.
- Do not wet the NOOĀNCE helmet, controller, or adapter. This is an electrical hazard that can result in personal injury. 
- Please do not use the NOOĀNCE helmet while charging the controller. This can lead to faults or electrical hazards. 

- Do not pull or forcibly untie the power cord when using the NOOÄNCE helmet. And this may cause a fault or electrical hazard.
- Please do not use it in the bathroom or the shower. Please do not use a NOOÄNCE helmet that has been wet with water.  
- Keep the NOOÄNCE helmet away from hot objects like ovens and indoor heaters.
- Keep your hair dry when using the device. Do not use the NOOÄNCE helmet when your hair is wet.
- Use only the controller and power adapter provided with the NOOÄNCE helmet.
- Use the NOOÄNCE helmet according to the instruction manual provided. Please do not use it for any other purpose.
- If you drop the NOOÄNCE helmet or suspect that the NOOÄNCE helmet is faulty (e.g., the LED or laser does not work, the NOOÄNCE helmet is cracked). Please get in touch with customer service for repair or replacement.
- Do not use the NOOÄNCE helmet when room temperature exceeds 33.3°C(92°F). This could cause dryness, redness or irritation of your scalp.
- If the cord is damaged on your NOOÄNCE helmet, do not use the NOOÄNCE helmet because you could get an electrical shock. Please get in touch with customer service for repair or replacement if the cord gets damaged.

- This NOOANCE helmet has no user-serviceable parts. Do not attempt to modify or repair this the NOOANCE helmet.
- Always unplug the NOOANCE helmet from the power outlet prior to cleaning the NOOANCE helmet. Failure to do so may result in an electrical hazard that may cause harm or injury.

Precaution

- If you are taking any medications or products that make you sensitive to light, you should test your NOOANCE helmet to make sure you do not have a reaction. Some of the things that can cause increased sensitivity to your NOOANCE helmet could be medicines for colds, allergies, pain medicine and drugs to treat infections.
- ⚠ Only use your NOOANCE helmet as directed. Do not try to change the way your NOOANCE helmet operates. Altering controls, adjustments, performance, or procedures other than those specified herein may result in exposure to hazardous radiation (laser). Never use acetone or other solvents on any part of your NOOANCE helmet. The use of solvents on the window can block the laser light from coming through.


Product Description

The wavelength of the NOOANCE LASER AND LED HELMET is a combination of 51 low-level lasers and 69 650nm light-emitting diodes.

Device Specifications

Product Name	LASER AND LED HELMET
Model No.	M-120
Output optical power density	<5mW
Wavelength	650±5nm
Number of Light Emitting Diodes (LED)	69
Number of Lasers (Class 3R)	51
Mode	Mode 1 is continuous mode; Mode 2 is pulse mode.
Operating temperature/humidity/pressure	41°F - 104°F (5°C -40°C), 15%-90% RH, 700 hpa to 1060 hpa
Storage & transportation temperature/humidity	-13°F -158°F (-25°C -70°C), 10%-90% RH
Ingress of water or particulate matter	IP22 (Enclose equipment without protection against ingress of water)
Charging supply	5V = 2A
Battery	3600mAh (Lithium-ion)
Treatment time	25 minutes
Charging time	4 hours
Recommended use	Three times a week on alternate days times
Dimensions of the device	26.2*21.7*14.5cm
Weight of the device	0.570kg / 20.11oz

Device Accessories

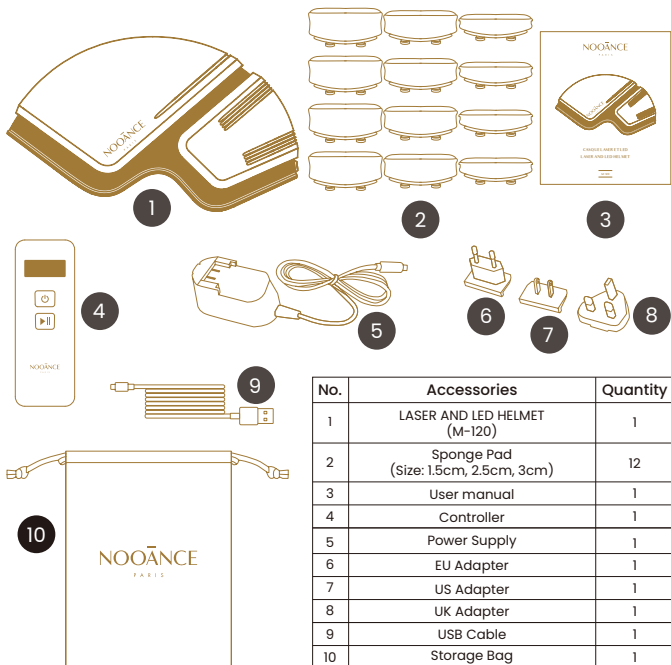
No.	Accessories	Quantity
1	LASER AND LED HELMET (M-120)	1
2	Sponge Pad (Size: 1.5cm, 2.5cm, 3cm)	12
3	User manual	1
4	Controller	1
5	Power Supply	1
6	EU Adapter	1
7	US Adapter	1
8	UK Adapter	1
9	USB Cable	1
10	Storage Bag	1

Light Sensitivity Test

We highly recommend users perform a light sensitivity test before their first use of the NOOÄNCE helmet.

1. Place the NOOÄNCE helmet on the inner forearm.
2. Turn on the NOOÄNCE helmet to start a treatment.
3. Wait until the treatment is done. The treatment will take 5 minutes to complete.
4. If your skin feels hot or uncomfortable during the treatment directly stop! You might be sensitive to light.
5. The NOOÄNCE helmet will automatically turn off after the treatment is done. Remove the NOOÄNCE helmet and wait for at least 24 hours.

Check your skin regularly over those 24 hours. If your skin appears red or itchy, the NOOÄNCE helmet is unsuitable for you. In the case of light sensitivity, please get in touch with customer support. If you do not see a reaction, proceed with the treatment instructions.

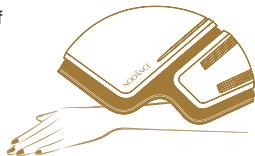
How To Use The Device


(Two ways to use)

General use

Preparation: Make sure your hair is dry before using the NOOANCE helmet.

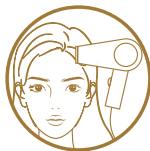① Connect the controller to use

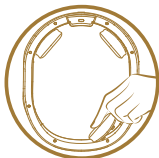Step 1

Please select the appropriate sponge pad size to fit inside the NOOANCE helmet.

Step 2

Connect the NOOANCE helmet to the controller using a USB cable.

Step 3

Position the NOOANCE helmet in its upright position on the top of your head.

Step 4

Sit down and position the dome comfortably. Long press the ON/OFF  Button to turn on the NOOANCE helmet.

(**Note:** If the start/pause button is not pressed within five minute, the NOOANCE helmet will turn off automatically.)

Step 5

Tap the ON/OFF  Button to select the preferred mode. (Mode 1 is continuity; Mode 2 is the pulse.)

Step 6

Press the start/pause button to begin the 25-minute treatment. A "beep" indicates the treatment has started. Once the treatment is complete, the NOOANCE helmet will turn off automatically with two "beep" sounds.

IMPORTANT: Unless the NOOANCE helmet is positioned on the user's head, it will not begin a treatment. If the NOOANCE helmet is removed from the head during a treatment, it will pause automatically. To resume the treatment, position the NOOANCE helmet on the user's head within 5 minutes and press the start/pause button.

Step 7

If you need to pause the treatment at any time, press the start/pause button, and the timer will stop. Press the start/pause button to resume the treatment. A start/pause audible "beep" sound will indicate the treatment has resumed. If the start/pause button is not pressed within 5 minutes during a paused treatment, the NOOANCE helmet will turn off automatically.

Step 8

To turn off the NOOANCE helmet, hold the button for 2 seconds.

IMPORTANT: The NOOANCE helmet has controller power outage memory. If an accidental USB cable unplugged occurs during treatment, the NOOANCE helmet will remember the treatment time. When you reinsert the USB cable into the controller, and the NOOANCE helmet is still located on the user's head, tap the start/pause button, and the NOOANCE helmet will automatically resume treatment. To turn the NOOANCE helmet off and start a new treatment at another time, hold the ON/OFF button for 2 seconds to turn off the NOOANCE helmet, then follow steps 2-8, listed above.

② Connect the adapter to use

Step 1

Connect the NOOANCE helmet to the adapter.

Step 2

Long press the On/Off button to turn on the NOOANCE helmet default start mode 1, and the indicator light on the NOOANCE helmet will always be white.

Step 3

Press the On/Off button again to switch to mode 2, and the indicator light on the NOOANCE helmet will flash white.

Step 4

Place the NOOANCE helmet vertically on your head after selecting your preferred mode. Once the treatment is complete, the NOOANCE helmet will turn off automatically with two "beep" sounds.

Step 5

To turn off the NOOANCE helmet, hold the button for 2 seconds.

Visibility of continuous mode and pulse mode

- Figure 1: Continuous mode. The LED and laser light up continuously, and the oscilloscope will display a fixed frequency.
- Figure 2: The pulse mode frequency is set to 1.97khz, so the rapid flicker of the laser is invisible to the naked eye. But it shows up clearly on the oscilloscope.

Figure 1: Mode 1 is continuous mode

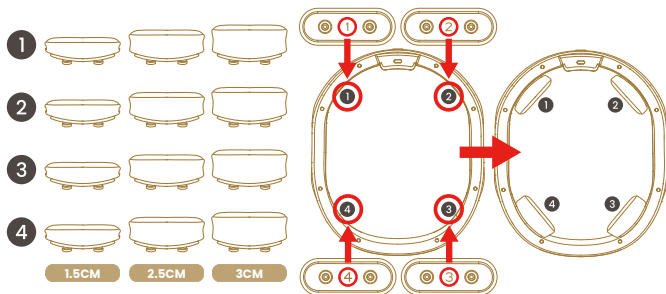
Figure 2: Mode 2 is pulse mode

Invisible pulse mode (2) is used in most laser machines or LEDs dedicated to professionals. It has the disadvantage of producing background noise associated with the pulse, which depends on each person's auditory sensitivity, which is why we have also created a continuous mode so that everyone can use our helmet.

Installation of bearings

There are 3 sizes to choose from. The corresponding number is located below each pad. Please put the correct number in the specified position.

Treatment Protocol

You should use the NOOANCE helmet for 25 minutes, every other day, on non-consecutive days.

Use the NOOANCE helmet for 16 weeks. Once the 16-week period is completed, use as frequently as needed to maintain the desired results without exceeding every other day for 25 minutes. Discontinuing the use of the NOOANCE helmet may result in the progression of hair loss.

Charging

Charge your NOOÄNCE helmet fully for 4 hours before your first use. Only use a 5V adapter or computer USB port (avoid fast chargers as they can damage the battery).

To charge the NOOÄNCE helmet insert the USB-C cable into the controller. Insert the USB-A side of the cable into a USB power supply. When the controller is charged, the LED indicator on the device flashes white, and when the controller is fully charged, the indicator will be steady white.

Before charging, make sure that the plug and socket are completely dry. Failure to do so may result in electric shock, short circuit, or fire.

Note: Always unplug the NOOÄNCE helmet before use, the NOOÄNCE helmet will not work while being charged.

Care And Storage

- Store the NOOANCE helmet in the supplied packaging. Make sure the NOOANCE helmet will not bend, stretch or crush while being stored or when travelled with.
- Keep the NOOANCE helmet away from sunlight and store the NOOANCE helmet in a dry location (e.g. away from windows, sinks, bathtubs, kitchens, and bathrooms) at around room temperature after each use.

Cleaning

- The NOOANCE helmet should be unplugged prior to any cleaning.
- Use a dry microfiber cleaning wipe or cloth to clean the device.
- For stains that cannot be removed with a dry cloth, a lightly water-moistened cloth may be used to clean the NOOANCE helmet.
- Do not use solvents to clean the NOOANCE helmet.

Troubleshooting

If the actions suggested fail to fix the problem, or the problem reoccurs, and your device is within the warranty period please return to NOOANCE.

⚠ DO NOT attempt to investigate any further.

⚠ DO NOT attempt to modify the NOOANCE helmet in any way. the NOOANCE helmet has NO parts you can repair yourself.

Problem	Cause	Solution
The NOOANCE helmet doesn't switch on.	The controller is not fully charged because the charging cable was not properly inserted into the power outlet.	Insert the charging cable properly into the power outlet.
	The controller is not fully charged because the charging cable was not properly inserted into the power outlet.	Insert the charging cable properly into the NOOANCE helmet and charge the NOOANCE helmet.
	The power outlet didn't work during charging.	Try a different 5V adapter or power outlet.
	The user does not wear the NOOANCE helmet upright on his/her head.	Unless the NOOANCE helmet is positioned on the user's head, it will not begin a treatment.
	The NOOANCE helmet connects not only to the controller but also to a power outlet.	After the controller is powered on, it is in charge mode by default. When the NOOANCE helmet cannot be activated in charging mode, disconnect the power, and the controller will resume regular use.

	Allow NOOANCE helmet internal software to reset itself.	Hold for 3-5 seconds and allow the NOOANCE helmet's internal software to reset itself.
The NOOANCE helmet is powered on but pressing does not activate the NOOANCE helmet.	The NOOANCE helmet will not power on unless the electronic sensor detects that the NOOANCE helmet is positioned on the users head. This is a built-in safety feature to protect the eyes from direct laser light exposure.	Activate the NOOANCE helmet by placing it vertically over your head.
The NOOANCE helmet is not lighting up, but the controller is on.	The power of the controller is low.	Please fully charge the controller before use.
	If the device still fails to light up after following the instructions and troubleshooting.	Please get in touch with customer service.
The NOOANCE helmet controller doesn't show the numbers clearly on the timer.	Try shutting down the device, restarting it, and allowing the NOOANCE helmet's internal software to reset itself.	If the timer cannot be displayed clearly after the NOOANCE helmet is restarted, please get in touch with customer service.

Electrical Safety

To reduce the risk of hazard by electric shock:

- Keep away from water. Contact with water may result in an electrical shock.
- Do not use while bathing.
- Do not place or store your NOOANCE helmet where it can fall or be pulled into a tub or sink.
- Do not expose your NOOANCE helmet to running water, or drop into water or any other liquid. If your NOOANCE helmet falls into water, unplug it immediately from the power outlet.
- If using your NOOANCE helmet near water, the use of an electrical outlet protected by a Ground Fault Circuit Interrupter is recommended.
- Do not dispose of your NOOANCE helmet in a fire or store in high temperatures greater than 140°F (60°C).
- Do not turn on or use your NOOANCE helmet while it is charging.
- The power supply cannot be separated during use. Do not position the equipment in a manner that will make it difficult to plug in or remove the power cord.
- The power cord is very long. To avoid danger, do not wrap the cord around your neck.

⚠ Basic safety precautions should be observed when using this, or any other electrical device, especially when children are present.

Travel and International Use

The NOOANCE helmet's lightweight design makes it a perfect travel companion. Here are some travel tips for you:

- If you are traveling to a different region from where you bought the NOOANCE helmet, please contact the retailer you purchased the helmet from that the NOOANCE helmet is only supplied with a 5V DC Low Voltage.
- Confirm with the steward that the NOOANCE helmet can be carried and/or used on the airplane prior to flight.

Disposal

DO NOT dispose of the device with household waste. The device should be recycled. The device contains a Lithium-ion battery. For protection of the environment, please safely dispose of the device following your local guidelines.

Warranty

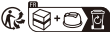
The NOOANCE helmet is warranted to be free from defects in material and workmanship for 24 months from its original purchase date under normal use. This warranty extends only to the original retail purchaser with original proof of purchase and only when purchased from an authorized retailer or reseller. This warranty only extends if the NOOANCE helmet is used in conjunction with authorized Device accessories. If the NOOANCE helmet should become defective within the warranty period, contact the retailer that you purchased the NOOANCE helmet from. This warranty does not cover the NOOANCE helmet damaged by:

- Accident, misuse, abuse, or alteration.
- Use with unauthorized accessories or use other than as instructed.
- Connecting to incorrect current and voltage.





Symbols

Electromagnetic compatibility complies with IEC 60601-1-2 Standard.

Symbole	Explication
	BF symbol, which indicates this device is in accordance to the degree of protecting against electric shock for type BF equipment.
	Indicates the Device is CE Marked and authorized to sell in Europe. (Lower-Voltage-Directive)
	The Waste of Electrical and Electronic Equipment (WEEE), Indicates the Device may not be discarded as household waste (2012/19/EU).
	Indicates that the electromagnetic radiation from the device is below the limits specified by the Federal Communications Commission and the manufacturer has followed the requirements of the Supplier's Declaration of Conformity authorization procedures.
	This warning symbol on the device label is a warning to read the enclosed documents before use and take all safety precautions required.
	Consult instructions for use.
	IP22
	Unique Device Identifier.
	This symbol means keep the Device dry.

	This symbol indicates the serial number of the device.
	This symbol refers to the batch code of the device.
	This symbol means keep the device away from direct sunlight.
	This symbol indicates the upper and lower operating temperatures that the device can be safely used at.
	This symbol means the upper and lower operating humidity that the device can be safely used at.
	This symbol means that the cardboard packaging that the device came in can be recycled.
	This symbol implies conformity to UK Directives on health, safety and environmental protection standards. UKCA stands for UK Conformity Assessed.
	Indicates the Device is CA65 Marked and authorized to be sold in the United States.
	Indicates the Device is ROHS Marked and authorized to sell in Europe.
	This symbol means that the device uses a Lithium-ion Polymer battery.
	The date of manufacture.

- Labels on package and device are in compliance with IEC 60825-1 Standard:

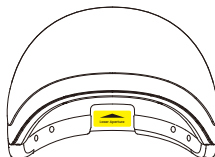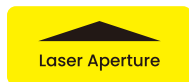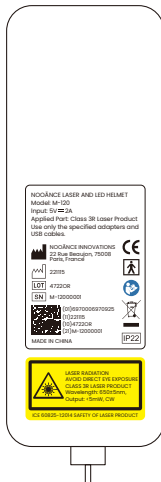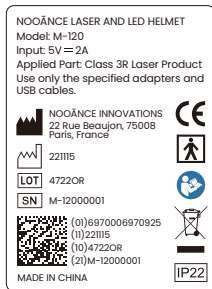NOOANCE INNOVATIONS
22 Rue Beaujon, 75008 Paris, France

DESIGNED IN FRANCE
MADE IN CHINA