
NOOĀNCE

P A R I S

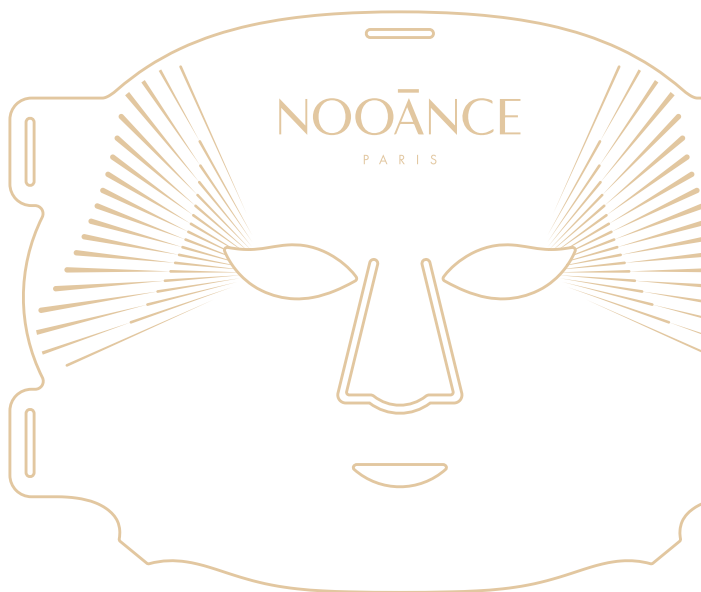MASQUE LED VISAGE

LED FACE MASK

PRO

Français ————— **03-16**

English ————— **18-31**

Sommaire

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Table of Contents

The science of light therapy	19
Contents of the box	20
Product Information	21
Device Specifications	
Frequency of use	
Charging the remote control/battery	
Light sensitivity test	23
Protective goggles	
Getting started	24
Installation of straps	
Remote control/Battery	
User Guide	
Battery operation	
Use of Mains Power	
Warnings	
Contraindications	
Care and storage	
Cleaning	
Warranty	
Air travel	
Disclaimer	
Symbols	30

The science of light therapy

The skin has the remarkable ability to absorb light energy. It uses this energy to push the cells to produce more energy (ATP), which increases their performance and allows them to carry out their natural functions, but at a faster rate.


Cells supercharged by light therapy use this energy to push the skin's fibroblasts to produce more collagen and elastin proteins and to stimulate the natural processes of cell renewal.

Your Nooance mask has 198 LEDs with an accuracy of +/- 2 nanometres.

- 99 red LEDs with a wavelength of 633 nanometres,
- 99 near-infrared LEDs with a wavelength of 830 nanometers. A synergy of 633+830nm for beneficial effects on the skin, stimulating various types of skin cells.

Contents of the box

1. NOOĀNCE LED Mask
2. Manual
3. Adjustable Straps
4. Power Supply
5. EU Adapter
6. US Adapter
7. UK Adapter
8. USB Cable
9. Silicone Eye Protectors
10. Silicone Eye Protectors
11. Remote Control/Battery
12. Storage Bag

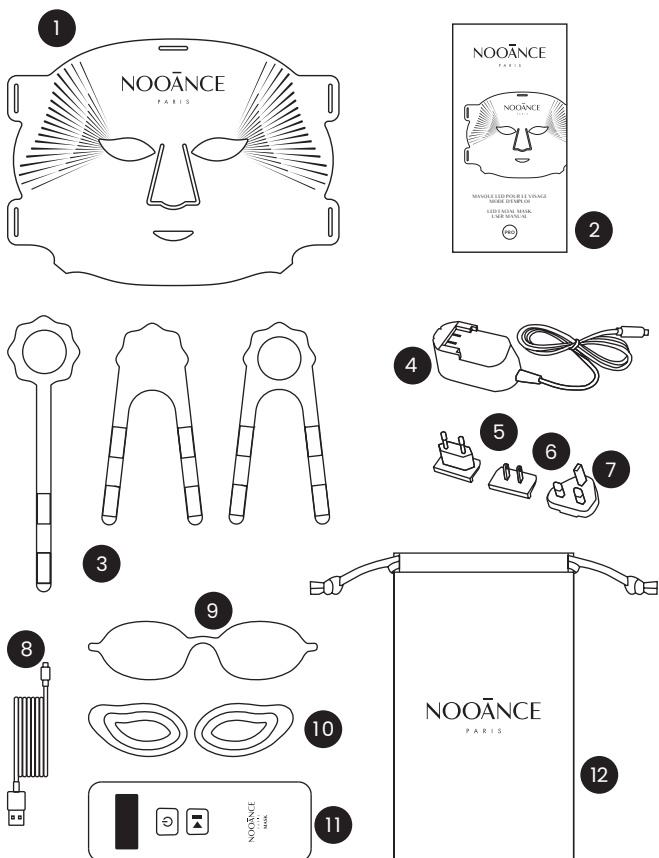Product Information

Device Specifications

LED wavelengths	633 nm +/- 2 nm (Red) 830 nm (Infrared)
Max. output optical power density	Red: 45mW/cm ² , Infrared: 45mW/cm ²
Energy delivered	2x27J/cm ²
Operating temperature	0°C to 40°C 32°F to 104°F
Storing and transport temperature	-5°C to 55°C 23°F to 131°F
Operating humidity	30% to 95% without condensation
Storage and Transport Humidity	20% to 95% without condensation
Charging supply	5V = 2A
Charging port	USB-C
Charging cable	USB-C to USB-A
Charging cable length	2 meters / 6.562 feet
Charging adapter power output	2A
Battery	3.7V, 3600mAh
Battery type	Lithium-ion
Charging time	240 minutes / 4 hours
Recommended use	30 to 60 minutes per week
Timer	10 minutes
Dimensions of the device	290*216*9 mm
Weight of the device	0.512 kg / 18.06 oz

The red light of the LEDs may appear orange after a few seconds of use as the eye adapts to this brightness. This is completely normal.

Frequency of use

We recommend using the mask for 30 minutes once a week (or 3 x 10 minutes spread over the week) for results on fine lines and wrinkles from as early as 16 weeks.

In the first month, you can start by doubling the time per week if you want to speed up the cell renewal process and boost skin firmness.

Charging the remote control/battery

Fully charge the remote control/battery for 4 hours before using for the first time.

Only use the provided equipment and avoid fast chargers as they may damage the battery.

Please use the adapter to connect the remote control for charging. Please note that the USB cable is only for connecting the LED mask, not for charging. The indicator light flashes during charging and becomes steady once the battery is fully charged.

Before charging the Remote Control/Battery, make sure that the plug and socket are completely dry. Failure to do so may result in an electric shock, short-circuit or fire.

NB: Unplug the remote control/battery from the mains before using it. The device does not work when it is charging.

Light sensitivity test

We strongly recommend that users carry out a light sensitivity test before first use.

1. Place the device on the inside of the forearm.
2. Switch on the device to start your session.
3. Turn the mask off after 5 minutes (the timer is set for 10 minutes by default).

If your skin becomes excessively hot or itchy during the session, stop immediately. You may be sensitive to light.

Check your skin regularly for the next 24 hours. If your skin is red and/or itchy, the appliance is not suitable for you. Please contact customer service.

Protective goggles

When using an LED mask, red light and infrared light pose no danger to the eyes. However, the use of goggles is recommended in the event of discomfort or hypersensitivity to light.

It is important not to use dark glasses (item 9) and silicone eye protectors (item 10) at the same time.

Using them together can cause uncomfortable pressure on the eyes. If the silicone protectors are adequate, they are sufficient; if the light remains too strong, it is preferable to opt for the dark glasses alone under the mask.

Getting started




Installation of straps


Sticky surface

Velvety surface

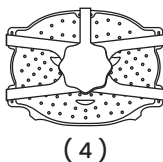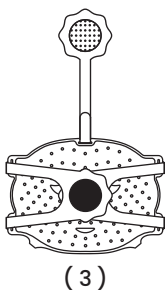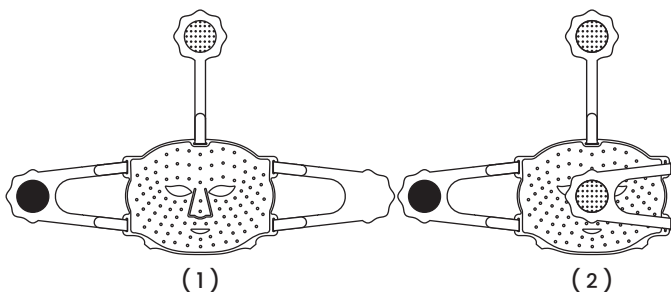Remote control/Battery

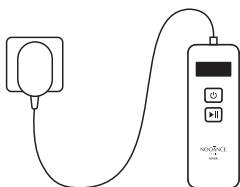
Charge your remote control/battery for 4 hours. You can then use it for 8 to 10 10-minute sessions.

User Guide

Preparation: Before starting the treatment, please thoroughly cleanse your skin.

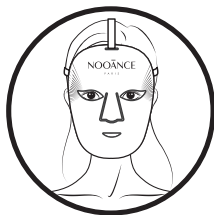
Apply the mask to your face, and strap it around your head using the adjustable straps.

To optimise the results of your Nooance mask, we recommend that you use one of our aqueous serums recommended on our website www.nooance-paris.com


At the end of the session, you can also apply a moisturising cream or serum

Battery operation

Connect the mask to the remote control/battery using the cable supplied (8)

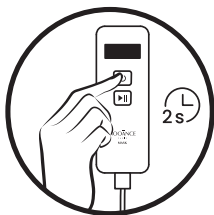Press and hold the ON/OFF button for 2 seconds to turn on the mask on the remote control.

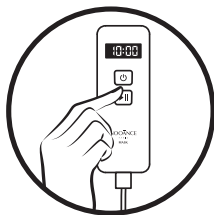Press the start/pause button to begin the 10-minute treatment. A "beep" indicates the treatment has started. Once the treatment is complete, the mask will turn off automatically with two "beep" sounds.

Use of Mains Power

Connect the mask directly to the mains using the 2m cable supplied and the mains adapter (you will need to time your session manually).

Do not connect the remote control/battery in this case.

⚠ Warnings

- Read all the instructions carefully before using this device.
- Stop using the device and consult your doctor if you experience any adverse reactions such as discomfort, dizziness, skin discolouration, skin irritation or any other reaction that seems abnormal.
- Stop using the device if it overheats or if you feel unpleasant heat on your skin.
- DO NOT allow children or people with reduced ability to play with the appliance. The appliance must not be cleaned or charged by children.
- DO NOT point the LED light directly into your eyes or those of another person.
- DO NOT use the appliance if it is not working properly, if it has been dropped or if it has been damaged. For warranty and repairs, contact customer service.
- DO NOT use, charge or store the appliance near hot surfaces, in contact with water or in direct sunlight.
- DO NOT attempt to repair the device. Repairing or opening the appliance will invalidate your warranty. The appliance contains no serviceable parts. This appliance contains non-replaceable batteries.
- DO NOT use the appliance with accessories not authorised by the manufacturer.
- DO NOT use solvents or abrasive cleaners on the appliance.
- DO NOT use the device in any context other than that specified in this manual.
- In very rare cases, temporary headaches may occur after a light therapy session. Consult your doctor if symptoms persist and discontinue use of the mask.
- Do not scratch or damage the light-emitting surface.
- Store the device as described in this manual.

- Inspect the device and do not use it if you see any signs of damage, such as a leaking battery. Contact our customer service department.
- Your warranty will be void if you do not use the device in accordance with the instructions in the manual.

Contraindications

- DO NOT use the device if you suffer from a photosensitive disorder (sensitivity to light) or if you are taking photosensitising medication or if you suffer from epilepsy. If in doubt, consult your doctor.
- DO NOT use the device if you have active cancer, a cancerous lesion, an open wound, warts or any skin condition. Consult your doctor if in doubt.
- DO NOT use the device on a bacterial infection, as this may aggravate the symptoms.
- DO NOT use the device if you are pregnant, planning to become pregnant or breastfeeding.
- DO NOT use the device on anyone under the age of 18.
- DO NOT use the device if you have a known allergic reaction to silicone. The device contains plastics and silicone that come into contact with the skin during use. If in doubt about allergies, contact your doctor before using the device.

Care and storage

- To maximise the life of the lithium battery, we recommend that you fully charge and discharge it during normal use.
- The device should be stored at room temperature in a dry place and in its original packaging to prevent it from deforming, stretching or being compressed. Avoid dusty, sunny or damp locations.

Cleaning

- Clean the device only with a clean, damp cloth or water wipe. You can also use an antibacterial cleaner or wipe.
- DO NOT clean the appliance when it is switched on or charging.

Warranty

- The LED mask comes with a 24-month guarantee applicable to the original purchaser and subject to proof of purchase, from an authorised dealer. This guarantee covers normal use with approved accessories. In the event of a fault, please do not hesitate to contact our customer service department.
- Please note that the guarantee does not cover damage caused by accidents, incorrect use, use of non-recommended accessories or poor electrical connection.

Air travel

The device's light weight makes it an ideal travel companion. It is allowed in the cabin. If you wish to check it into the hold, take the Remote Control/Battery with you in the cabin. The main airlines refuse to accept lithium batteries in the hold. if in doubt, do not hesitate to consult your airline.

Disclaimer

The NOOANCE Face Mask is not intended for use in the diagnosis of disease and other conditions, nor for the cure, mitigation, treatment or prevention of disease, nor is it intended to alter the structure or function of the body.

Troubleshooting

If the actions suggested fail to fix the problem, or the problem reoccurs, and your device is within the warranty period please contact our customer service by writing to bonjour@nooance-paris.com.

⚠ DO NOT attempt to investigate any further.

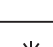
⚠ DO NOT attempt to modify the device in any way. The device has NO parts you can repair yourself.

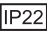
Problem	Cause	Solution
The device doesn't switch on.	The device is not fully charged, because the charging cable was not properly inserted into the power outlet.	Insert the charging cable properly into the power outlet.
	The device is not fully charged, because the charging cable was not properly inserted into the device.	Insert the charging cable properly into the device and charge the device.
	The power outlet didn't work during charging.	Try a different 5V adapter or power outlet.
	The on/off button has not been pushed 2 seconds full second.	Push the on/off button for 2 seconds straight.
The device switches off unexpectedly.	The battery level is too low.	Charge the Remote Control/Battery.

⚠ Disposal

DO NOT dispose of the device with household waste. The device should be recycled. The device contains a Lithium-ion battery. For protection of the environment, please safely dispose of the device following your local guidelines.

Symbols

Symbol	Explanation
	Indicates the device is CE Marked and authorized to sell in Europe. (Lower-Voltage-Directive)
	The Waste of Electrical and Electronic Equipment (WEEE), Indicates the device may not be discarded as household waste (2012/19/EU).
	Indicates that the electromagnetic radiation from the device is below the limits specified by the Federal Communications Commission and the manufacturer has followed the requirements of the Supplier's Declaration of Conformity authorization procedures.
	This warning symbol on the device label is a warning to read the enclosed documents before use and take all safety precautions required.
	Consult instructions for use.
	Date of Manufacture.
	This symbol means keep the device dry.
	This symbol indicates the serial number of the device.
	This symbol refers to the batch code of the device.
	This symbol indicates the USB connector.
	This symbol means that the device uses a Lithium-ion Polymer battery.
	This symbol means keep the device away from direct sunlight.
	This symbol means the upper and lower operating temperatures that the device can be safely used at.
	This symbol indicates the upper and lower operating humidity that the device can be safely used at.

	<p>This symbol means that the cardboard packaging that the device came in can be recycled.</p>
	<p>IP22</p>
	<p>Type BF Applied Part, B indicates it is used for the body and F indicates that it is not connected to the mains voltage.</p>
	<p>Unique Device Identifier.</p>
	<p>This symbol implies conformity to UK Directives on health, safety and environmental protection standards. UKCA stands for UK Conformity Assessed.</p>
	<p>Indicates the device is CA65 Marked and authorized to sell in the United States.</p>
	<p>Indicates the device is ROHS Marked and authorized to sell in Europe.</p>

Designed in France

Made in China

NOOANCE INNOVATIONS
 22 Rue Beaujon, 75008 Paris, France
bonjour@nooance-paris.com